


pTurboRFP-mito vector

The vector sequence has been compiled using the informa-tion from sequence databases, published literature, and other sources, together with partial sequences obtained by Evrogen. This vector has not been completely sequenced.

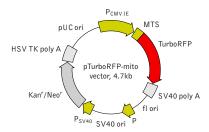

For vector sequence, please visit our Web site at http://www.evrogen.com/products/vectors.shtr

P_{CMV IE}: 1-589 Enhancer region: 59-465 TATA box: 554-560 Transcription start point: 583

Lo

TurboRFP-mito fusion Start codon (ATG): 597-599

Mitochondrial targeting sequence (MTS): 597-683 Start of TurboRFP coding sequence (ATG): 705-707 Stop codon (TGA): 1398

SV40 early mRNA polyadenylation signal Polyadenylation signals: 1554-1559 & 1583-1588

mRNA 3' ends: 1592 & 1604 f1 single-strand DNA origin: 1651-2106

Bacterial promoter for expression of Kan^r gene -35 region: 2168-2173; -10 region: 2191-2196 Transcription start point: 2203

SV40 origin of replication: 2447-2582

SV40 early promoter Enhancer (72-bp tandem repeats): 2280-2351 & 2352

21-bp repeats: 2427-2447, 2448-2468 & 2470-2490

Early promoter element: 2503-2509 Major transcription start points: 2499, 2537, 2543 & 2548

Kanamycin/neomycin resistance gene Neomycin phosphotransferase coding sequences: Start codon (ATG): 2631-2633; Stop codon: 3423-3425 G->A mutation to remove Pst I site: 2813

C->A (Arg to Ser) mutation to remove BssH II site: 3159 Herpes simplex virus (HSV) thymidine kinase (TK) polyadenylation signal

Polyadenylation signals: 3661-3666 & 3674-3679 pUC plasmid replication origin: 4010-4653

Product	Cat.#	Size	
pTurboRFP-mito vector	FP237	$20~\mu \mathrm{g}$	
Vector type	mammalian ovn	rossion vootor	
· ·	mammalian expression vector		
Reporter	TurboRFP		
Reporter codon usage	mammalian		
Promoter for TurboRFP	P _{CMV IE}		
Host cells	mammalian		
Selection	prokaryotic - kanamycin		
	eukaryotic - neomycin (G418)		
Replication	prokaryotic - pUC ori		
•	eukaryotic - SV40 ori		
Use	red (orange) fluorescent labeling of mitochondria		

scription

pTurboRFP-mito is a mammalian expression vector intended for red (orange) fluorescent labeling of mitochondria in living cells. The vector encodes red (orange) fluorescent protein TurboRFP fused to mitochondrial targeting sequence (MTS) derived from the subunit VIII of human cytochrome C oxidase [Rizzuto et al. 1989; Rizzuto et al. 1995]. MTS is fused to the TurboRFP N-terminus.

TurboRFP codon usage is optimized for high expression in mammalian cells (humanized) [Haas et al. 1996].

pTurboRFP-mito vector can be used as a source of TurboRFP-MTS hybrid sequence. The vector backbone contains unique restriction sites that permit its excision and further insertion into expression vector of choice.

Note: The plasmid DNA was isolated from dam +-methylated E.coli. Therefore some restriction sites are blocked by methylation. If you wish to digest the vector using such sites you will need to transform the vector into a dam host and make fresh DNA.

The vector backbone contains immediate early promoter of cytomegalovirus (P_{CMV IE}) for protein expression, SV40 origin for replication in mammalian cells expressing SV40 T-antigen, pUC origin of replication for propagation in E. coli, and f1 origin for single-stranded DNA production. SV40 polyadenylation signals (SV40 poly A) direct proper processing of the 3'-end of the reporter mRNA.

SV40 early promoter (PSV40) provides neomycin resistance gene (Neor) expression to select stably transfected eukaryotic cells using G418. Bacterial promoter (P) provides kanamycin resistance gene expression (Kan^r) in E. coli. Kan^r/Neo^r gene is linked with herpes simplex virus (HSV) thymidine kinase (TK) polyadenylation signals.

Expression in mammalian cells

pTurboRFP-mito vector can be transfected into mammalian cells by any known transfection method. CMV promoter provides strong, constitutive expression of the TurboRFP-MTS fusion in eukaryotic cells. If required, stable transformants can be selected using G418 [Gorman 1985].

Propagation in E. coli

Suitable host strains for propagation in E. coli include DH5alpha, HB101, XL1-Blue, and other general purpose strains. Plasmid incompatibility group is pMB1/CoIE1. The vector confers resistance to kanamycin (30 μ g/ml) to E. coli hosts. Copy number in E. coli is about 500.

References

Gorman, C. (1985). "High efficiency gene transfer into mammalian cells." In: DNA cloning: A Practical Approach, Vol. II. Ed. by Glover. (IRL Press, Oxford, U.K.) Pp. 143-190.

Haas, J. et al. (1996) "Codon usage limitation in the expression of HIV-1 envelope glycoprotein." Curr Biol, 6 (3): 315-324 /

Rizzuto, R. et al. (1989) "A gene specifying subunit VIII of human cytochrome c oxidase is localized to chromosome 11 and is expressed in both muscle and non-muscle tissues." J Biol Chem, 264 (18): 10595-10600 / pmid: 2543673

Rizzuto, R. et al. (1995) "Chimeric green fluorescent protein as a tool for visualizing subcellular organelles in living cells." Curr Biol, 5 (6): 635-642 / pmid: 7552174

Notice to Purchaser:

TurboRFP-related materials (also referred to as "Products") are intended for research use only.

The Products are covered by U.S. Pat. 8,138,320; European Pat. 1994149; and other Evrogen Patents and/or Patent applications pending. By use of these Products, you accept the terms and conditions of the applicable Limited Use Label License #001: http://www.evrogen.com/products/Evrogen-FP-license.shtml.

The CMV promoter is covered under U.S. Patents 5,168,062 and 5,385,839, and its use is permitted for research purposes only. Any other use of the CMV promoter requires a license from the University of Iowa Research Foundation, 214 Technology Innovation Center, Iowa City, IA 52242